




I-TECH Lite



Brushless Micromotor

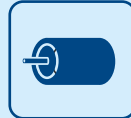

**Autoclavable Brushless Micromotor
Guaranteed up to 1000 cycles**



Electronic Calibration of ContraAngle



40000 RPM



Torque: 72 Ncm



Touch Screen



**CA
20:1**



**Available with
Conrangle 20:1 Fiber Optic**

LED

**High Luminesity
Light**




I-TECH Lite

Implant Motor with **4.3" TouchScreen/LED Lite**
brushless motor - autoclavable
Motor guaranteed 1000 sterilization cycles
Form and Functions perfectly mixed



I-TECH Lite

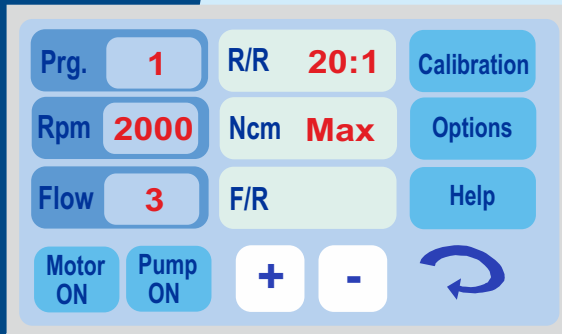

LED

- Wide light area
- Light without shadows
- Power: 31000Lux

I-TECH LITE is a complete system, specially developed under consideration of all the needs in implantology and surgery: **design, TouchScreen LED LITE, powerful performance and precise functionality.**

Display TouchScreen 4.3"

4.3" glass touch-screen intuitive control panel for user-friendly piloting of all programs. Select personalised settings, display information just the way you want and enjoy simple yet precise implantology unit control

Powerful Motor - L.E.D. Lighting

72 Ncm (at the instrument) of pure power for all clinical performances

STERILIZABLE: 135°C (275°F) **guaranteed up to 1000 sterilisation cycles.**

ATS LED's generate natural daylight-quality light to illuminate the treatment area enabling more precise surgery and shortened operation times. **Constant intensity at both High and Low speed - White Light - Clear color contrasts.**


Select the Programs

«DeLuxe» Multifunction Foot Control

Set the peristaltic pump flow

Select FWD/REV

Turn ON Motor/Pump

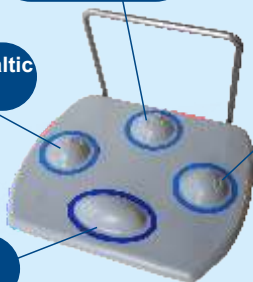


Covered entirely with syntetic rubber: resistant to any leaking fluid. The reostatic effect is ensured by means of a **specialized «pressure sensor»**. There are NO mechanical or moving parts.



Advanced Irrigation Pump

- Easy set up of irrigation tubes - Very ergonomic opening and closing system - Better control of irrigation flow with very constant flow rate - Stable and quiet operations - 3 Flow Rates available: MIN - MED - MAX

ATS Dental

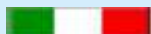

Via Vecchia Provinciale Lucchese, 49
51034 SERRAVALLE PISTOIESE (PT)-ITALY
info@atsdental.it - www.atsdental.it

Model	Code	Description
I-TECH LITE	ITL500	Implantology Unit: - Electronic Control Console - Autoclavable Brushless Micromotor LED Lite - «DeLuxe» Foot Control - Solution Holder - Motor Holder - 1 Tube