

«Omega» Implant «Plus»



Brushless Micromotor



Autoclavable Brushless Micromotor
Guaranteed up to 1000 cycles



Electronic Calibration of ContraAngle



40000 RPM




Torque: 72 Ncm



The new generation...
of Power and Control

**CA
20:1**

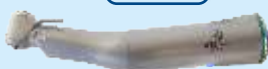

Available with
Constrangle 20:1 Fiber Optic

LED

High Luminosity
Light




«Omega» Implant «Plus»

Our powerful Implant Motor with integrated coolant pump is easy and simple to use - **LED light**
Reliable and convenient **Autoclavable Brushless motor**



LED
- Wide light area
- Light without shadows
- Power: 31000Lux

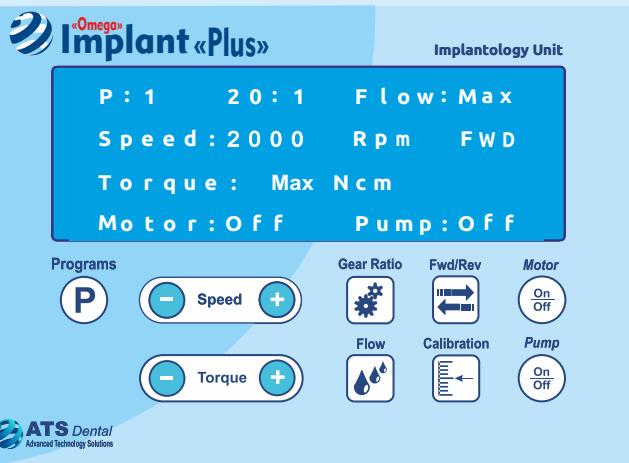


Our **OMEGA IMPLANT «Plus»** is a complete system, specifically designed with all the power and functionality needed for intelligent and precise implantology and surgery: **powerful performance and precise functionality**. See for yourself!

Complete, powerful and precise!

Our **OMEGA IMPLANT «Plus»** is a compact tabletop unit for implantology that combines the power of a **BRUSHLESS motor** with a user friendly . All necessary treatment stages of implantology can be perormed safely and precisely **using just one 20:1 contra-angle**.

OMEGA IMPLANT «Plus» is equipped with **Contra-Angle calibration**. The motor torque is **72 Ncm** and the motor speed range is **400 – 40,000 rpm (20 - 2000 with 20:1 contra-angle)** .

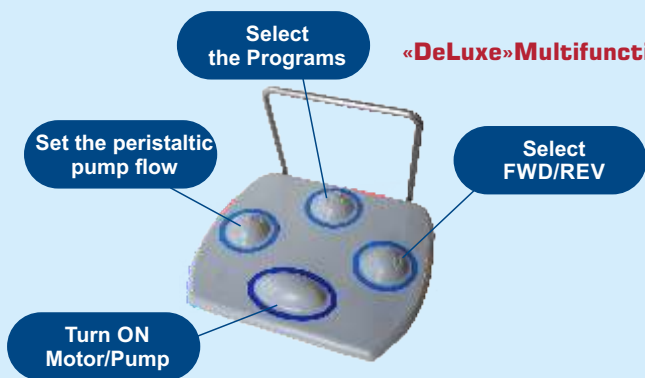

«DeLuxe» Multifunction Foot Control



Advanced Irrigation Pump

Covered entirely with syntetic rubber: resistant to any leaking fluid. The reostatic effect is ensured by means of a **specialized «pressure sensor»**. There are **NO** mechanical or moving parts.

Easy set up of irrigation tubes - Very ergonomic opening and closing system - Better control of irrigation flow with very constant flow rate - Stable and quiet operations - **3 Flow Rates available: MIN - MED - MAX**



Via Vecchia Provinciale Luccese, 49
51034 SERRAVALLE PISTOIESE (PT)-ITALY
info@atsdental.it - www.atsdental.it

Model	Code	Description
OMEGA IMPLANT Plus	AOMG300DL	Implantology Unit: - Electronic Control Console - Autoclavable Brushless Micromotor - «DeLuxe» Foot Control - Solution Holder - Motor Holder - 2 Tube