

TAC UltraMix



fully computerized
and programmable capsule mixer



TAC UltraMix

SIMPLY TO USE

The touch control system is quick and simple to operate. The capsule holder design makes loading and removal easy

TOUCH POAD SELECTOR

4200 rpm speed range ensures the most complete mixing possible

UNIVERSAL ARM DESIGN

Multi use mixer suitable for all types of encapsulated dental materials

REINFORCED CARBON-FIBER ARM

More durable than metal arms

VIBRATION ABSORBERS.

QUITE OPERATION AND LOW VIBRATION
Minimal vibration of all external parts ensures it will not move from its standing position

PROTECTIVE COVER


Covers arms and capsule during high speed mixing to reduce noise. It will immediately stop mixing when the protective cover is open

RAMP TO SEEK QUICKLY

It goes to 0 to set cycles per minute immediately

EASY TO DISINFECT

The hygienic one surface covers facilitate infection control cleaning and the non-porous plastic surface can be cleaned with any dental disinfectant.

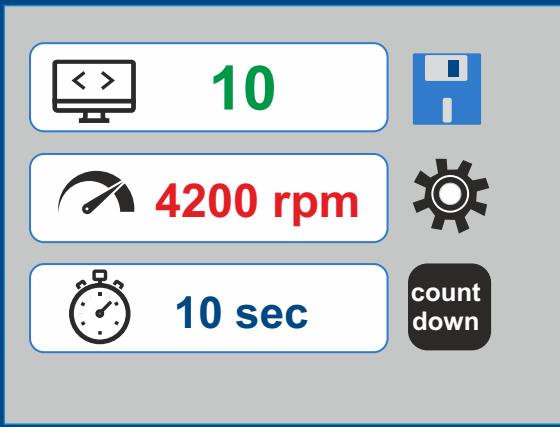

TAC UltraMix is a modern universal mixing device for amalgam, glass ionomer cements and other predosed dental materials in capsules and injection capsules.

This well-engineered and powerfull unit has a reliable «speed and time control».

The microprocessor controls:

- «10 Programs»
- «Speed»: from 3000 rpm to 6000 rpm
- «Time»: from 0 to 99 seconds

Mixing.Strong.Always

fully computerized
and programmable
capsule mixer

Model	TAC-UltraMix	Class MD	I rule 12 in Annex IX 93/42/EEC
Dimensions	300x220x170	Console	IP21
Materials	ABS	7488 Classification	Type 2
Noise	< 65 dBA	Operating Mode	30" ON – 180" OFF
Supply Voltage	AC 100-240V out 30V	Max Working Time	30'
Frequency	50/60 Hz	Fuses	NA

ATS Dental



Via Vecchia Provinciale Lucchese, 49
51034SERRAVALLE PISTOIESE (PT)
ITALY
info@atsdental.it - www.atsdental.it

