

TAC UltraMix

programmable capsule mixer


LINEA TAC

AMALGAMATOR

Mixing.Strong.Always

Features and Benefits

TAC UltraMix is a new well-engineered amalgamator. This powerful unit has a reliable «**speed and time control**». The microprocessor controls «**10 Programs**», the «**Speed**» - that can be adjusted from 2000 rpm to 7000 rpm- and the «**Time**» - that can be adjusted from 0 to 99 seconds. Cycles per minute are continuously verified.

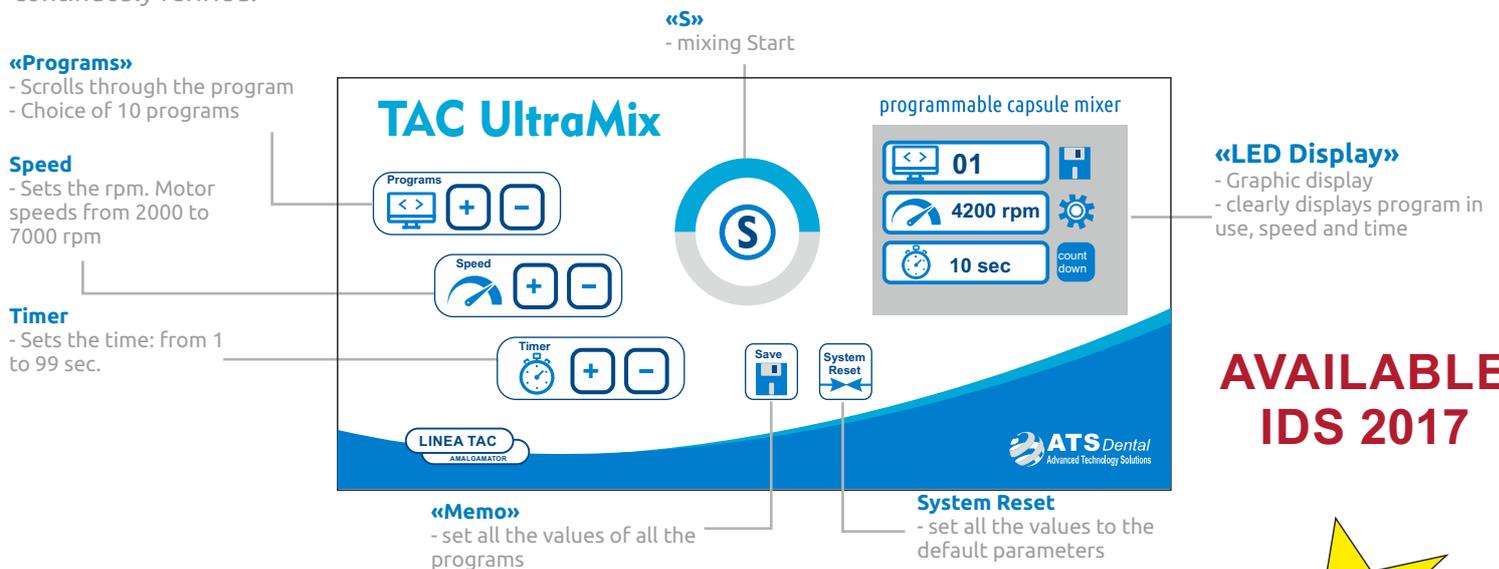

AVAILABLE
IDS 2017



Simply to use

The touch control system is quick and simple to operate. The capsule holder design makes loading and removal easy

Touch Pad Selector

Allow easy time setting

Universal Arm Design

Multi use mixer suitable for all types of encapsulated dental materials

Reinforced Carbon-Fiber Arms

More durable than metal arms

Vibration absorbers. Quite operation and low vibration

Minimal vibration of all external parts ensures it will not move from its standing position

Protective Cover

Covers arms and capsule during high speed mixing to reduce noise. It will immediately stop mixing when the protective cover is open

Ramp to seed quickly

Goes to 0 to 4200 cycles per minute immediately

Easy to disinfect

The hygienic one surface covers facilitate infection control cleaning and the non-porous plastic surface can be cleaned with any dental disinfectant.



TECHNICAL DATA

Model	TAC UltraMix	Class MD	I rule 12 in Annex IX 93/42/EEC
Dimensions	300x220x170	Console	IP21
Materials	Metal	7488 Classification	Type 1
Noise	< 65 dBA	Operating Mode	30" ON - 180" OFF
Supply Voltage	AC 230V	Max Working Time	30'
Frequency	50/60 Hz	Fuses	1xT2.0A 250V

CE  Made in Italy

ATS Dental
Advanced Technology Solutions